


SafeSheath® II

HEMOSTATIC PEEL AWAY INTRODUCER SYSTEM
FOR VASCULAR ACCESS

LATEX-FREE

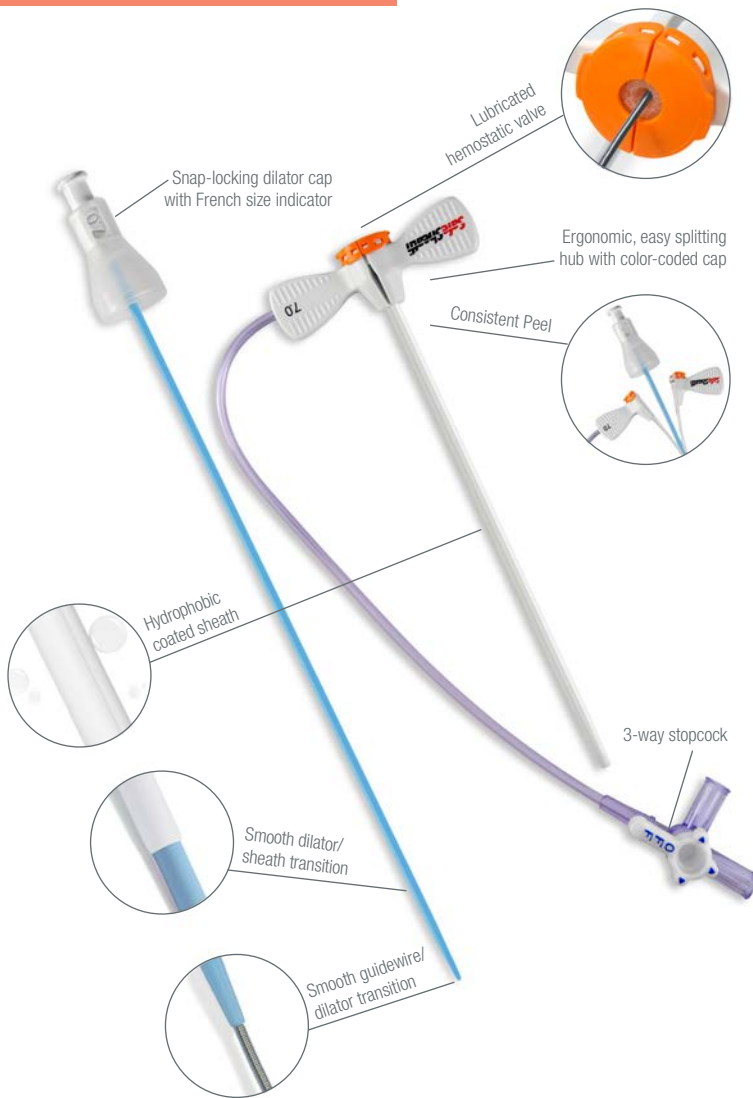

Introducing the **Adelante® SafeSheath® II Hemostatic Peel Away Introducer System for Vascular Access**. This next generation peel away introducer features the latest lubricated hemostatic valve membrane that provides low insertion forces during procedures. It also includes a sideport with a 3-way stopcock that provides a convenient means of aspirating and flushing the introducer.

Adelante® SafeSheath® II with Side Port Features:


- Lubricated hemostatic valve to reduce the risk of blood loss and air embolism
- Sideport with 3-way stopcock for infusion and contrast injection
- Consistently clean and smooth peel
- Ergonomic, easy-splitting break-away hub
- French size printed on handle and indicated by color-coded introducer cap
- Snap-locking dilator connector
- Accepts maximum guidewire up to 0.038"/0.96 mm

Model # Kit	French Size	Sheath Length	Dilator Length
<i>(Five per box)</i>			
STANDARD			
AS06013B5	6 F	13 cm	19.5 cm
AS07013B5	7 F	13 cm	19.5 cm
AS08013B5	8 F	13 cm	19.5 cm
AS85013B5	8.5 F	13 cm	19.5 cm
AS09013B5	9 F	13 cm	19.5 cm
AS95013B5	9.5 F	13 cm	19.5 cm
AS10013B5	10 F	13 cm	19.5 cm
AS10513B5	10.5 F	13 cm	19.5 cm
AS11013B5	11 F	13 cm	19.5 cm
AS12013B5	12 F	13 cm	19.5 cm
AS12513B5	12.5 F	13 cm	19.5 cm
LONG			
AS06023B5	6 F	23 cm	29.5 cm
AS07023B5	7 F	23 cm	29.5 cm
AS08023B5	8 F	23 cm	29.5 cm
AS09023B5	9 F	23 cm	29.5 cm
AS10023B5	10 F	23 cm	29.5 cm
AS10523B5	10.5 F	23 cm	29.5 cm
AS11023B5	11 F	23 cm	29.5 cm

Adelante® SafeSheath® II Kit Contains:

- One SafeSheath® Ultra peel-away sheath
- One dilator
- One guidewire with a flexible straight and J tip with straightener
- One needle with protector
- One syringe

*Introducers supplied sterile. For single use only.
Standard lengths are noted above. Other lengths available upon request.
For additional information e-mail us at sales@oscor.com*

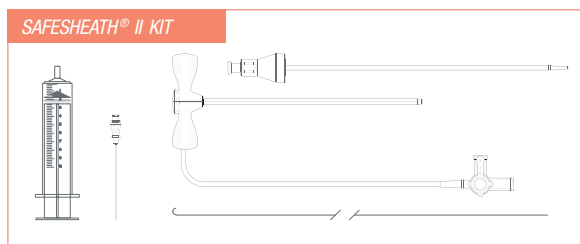

H. Lee SafeSheath

** SafeSheath is a registered trademark of Pressure Products, Inc.*



P: 727.937.2511
F: 727.934.9835
E: sales@oscor.com
www.oscor.com

CAUTION: Federal Law restricts this device to sale by or on the order of a physician (U.S.A). Please refer to the Instructions for Use for a complete listing of the indications, contraindications, precautions and warnings, where applicable.

