

Portal Vein Access

- **Scorpion® Portal Vein Access Set**
- **Scorpion®X Portal Vein Access Set**
- **Traveler™ Portal Vein Access Sets**
- **Intara™ Introducer Sheath**

Scorpion® and Scorpion®X Portal Vein Access Set

The Scorpion® Portal Vein Access Series is designed to redefine procedural control through dynamic steering and enhanced directional precision.

Greater Procedure Visibility


- Echogenic tips enable visibility of the needle or stylet on various imaging modalities for easy detection while piercing through liver to portal vein
- Alignment marker helps identify when the needle/stylet is flush with the end of the stiffening cannula
- Radiopaque bands on the stiffening cannula provide amplified visibility under fluoroscopy to assist with locating and orienting each component

Improved Durability

Needle/Stylet catheters are constructed of PEEK material instead of traditional PTFE to reduce the risk of fishmouthing and buckling when advancing through cirrhotic parenchyma.

Unparalleled Steerability


- Independently steer the curved needle/stylet and the curved stiffening cannula to expertly maneuver through the vasculature and parenchyma and to change trajectory if necessary
- Directional handles allow precise directional control of the stiffening cannula and needle/stylet to enter the portal vein and change direction when advancing thorough the parenchyma

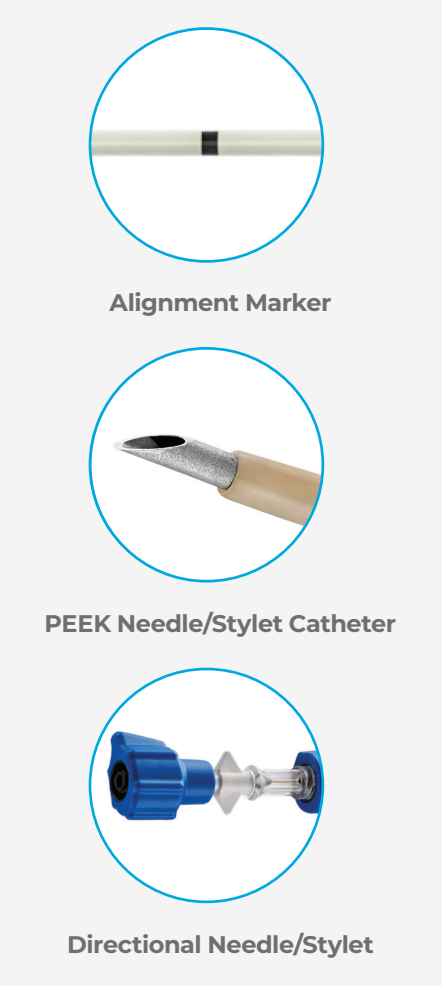

Box of 1		Cannula		Sheath		Needle		Needle Catheter		MPA Catheter	
Catalog #	Description	Gauge	Length	Gauge	Length	Diameter	Length	Gauge	Length	Gauge	Length
TPS005	Scorpion	14G	52cm	7F	52cm	.040"	73cm	5F	71cm	5F	80cm
TPS006	Scorpion X	13G	52cm	8F	52cm	17G	73cm	6.2F	71cm	5G	80cm


Each set includes a spacer clip

Scorpion and Scorpion X Portal Vein Access Set

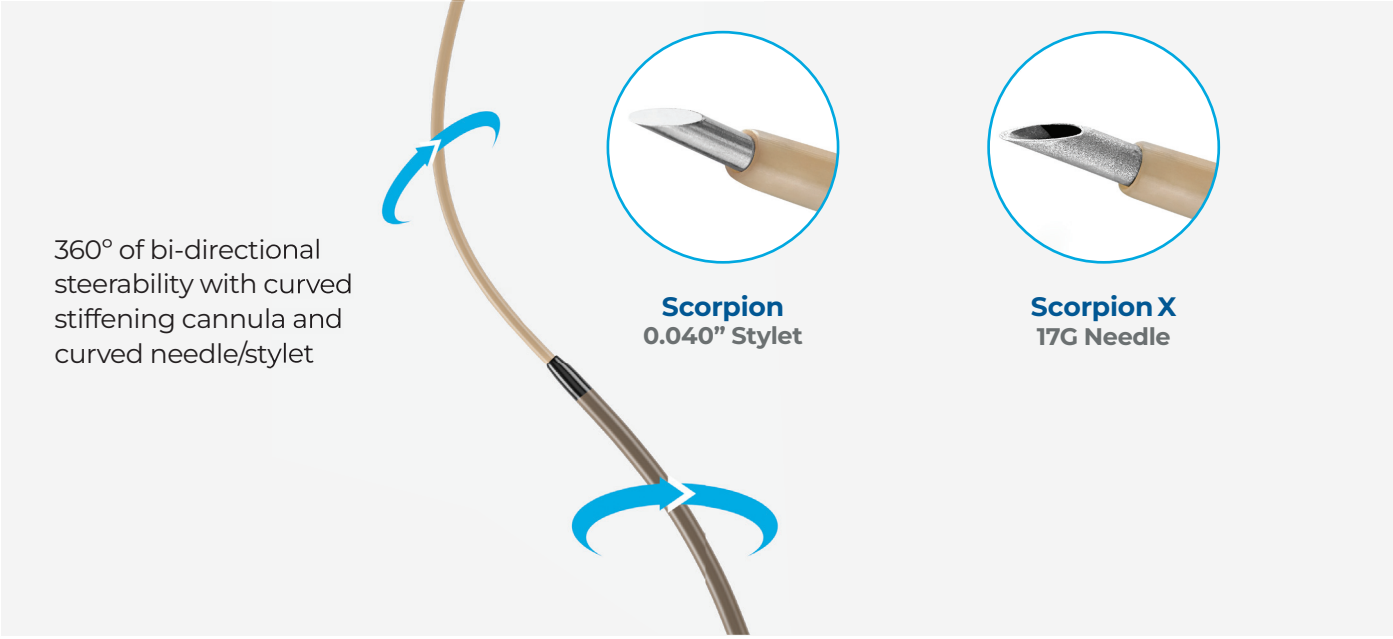

Traveler™ Portal Vein Access Sets

The Traveler™ Portal Vein Access Series provides superior visibility and durability, while maintaining current procedural techniques.

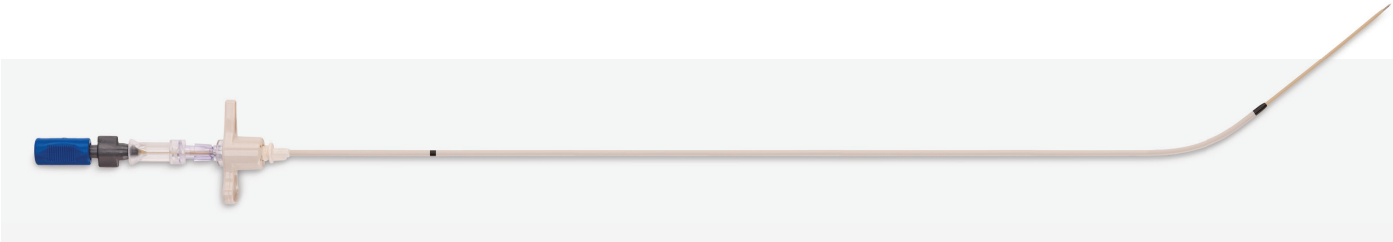

Box of 1		Cannula		Sheath		Needle		Needle Catheter		MPA Catheter	
Catalog #	Description	Gauge	Length	Gauge	Length	Diameter	Length	Gauge	Length	Gauge	Length
TPS001	Traveler 38	14G	52cm	7F	52cm	.038"	69cm	5F	66cm	5F	80cm
TPS002	Traveler 21	14G	52cm	7F	52cm	21G	72cm	5F	67cm	5F	80cm
TPS003	Traveler 16					16G	57cm	7F	51cm	5F	80cm

Each set includes a bending tool and spacer clip